

Resistivity Matrix Card Configuration – Keithley 7709


Wiring Color Code

- A Red
- B Black
- C Blue
- D Green
- I White
- V Yellow

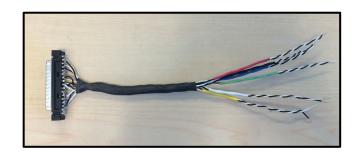


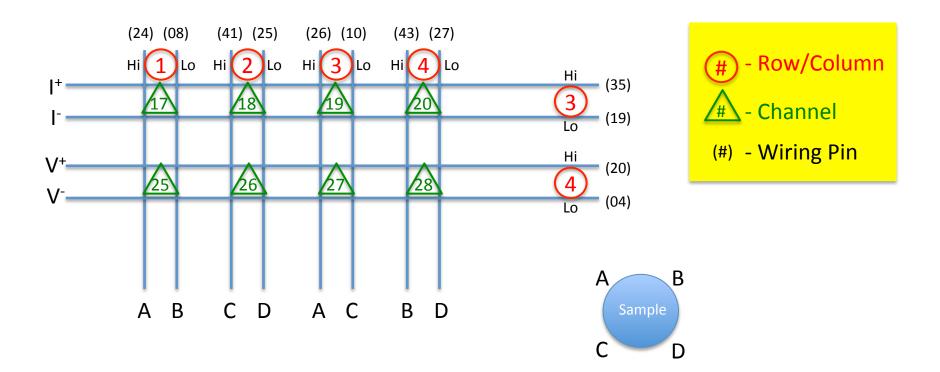


Steps 1, 2 Measurement Sequence – R_{AB, CD}

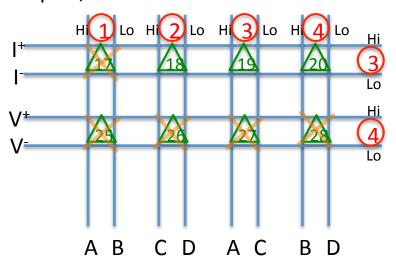

- Current Off
- 2. Short Card Close 17, 25, 26, 27, 28
- 3. Prepare Measurement Open 25, 27, 28
- 4. Current On, Delay, Measure Voltage
- 5. Current Reverse, Delay, Measure Voltage

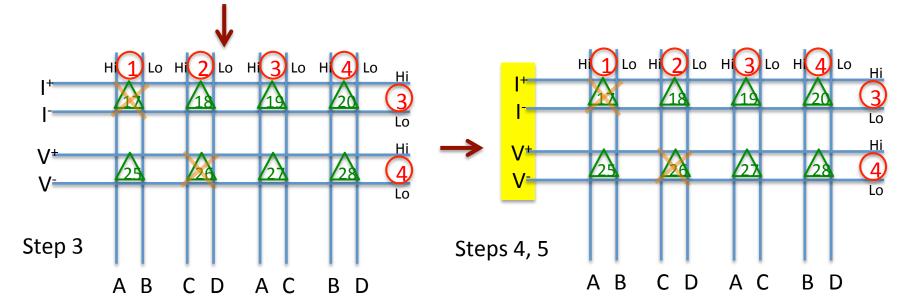


Measurement Sequence – R_{AB, CD}

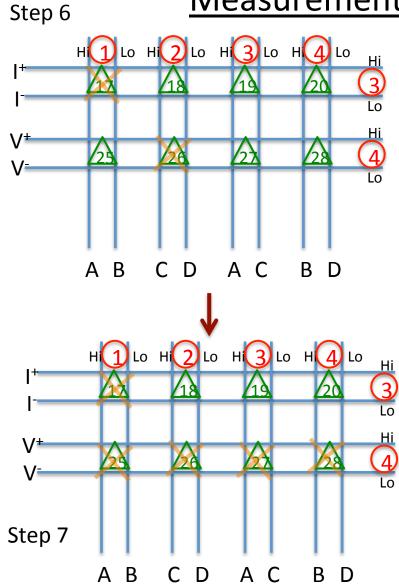


- 6. Current Off
- 7. Short Card Close 25, 27, 28


Steps 1, 2 Measurement Sequence – R_{CD, AB}

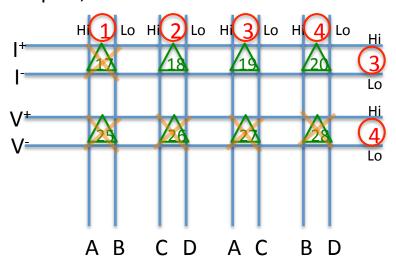


- Current Off
- 2. Short Card Close 17, 25, 26, 27, 28
- 3. Prepare Measurement Close 18
- 4. Prepare Measurement Open 17, 26, 27, 28

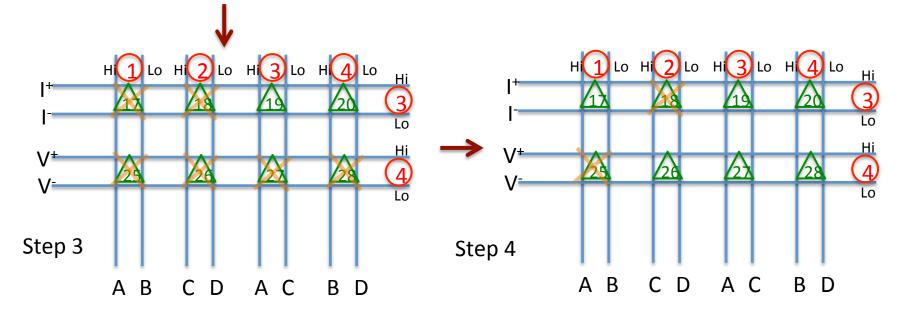


Steps 5, 6, 7 Measurement Sequence – R_{CD, AB}

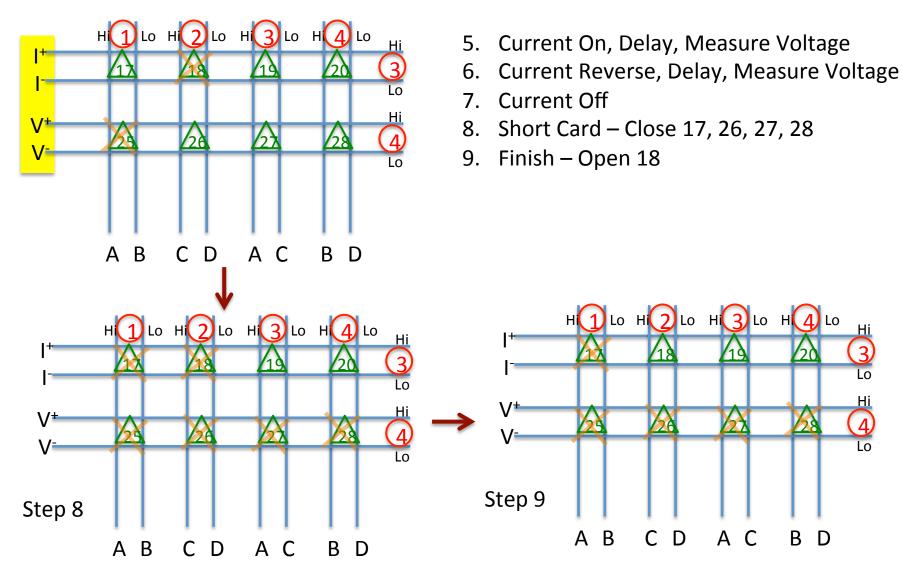


Steps 1, 2 Measurement Sequence – R_{AC, BD}

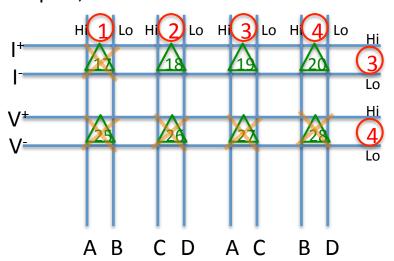


- Current Off
- 2. Short Card Close 17, 25, 26, 27, 28
- 3. Prepare Measurement Close 19
- 4. Prepare Measurement Open 17, 25, 26, 27

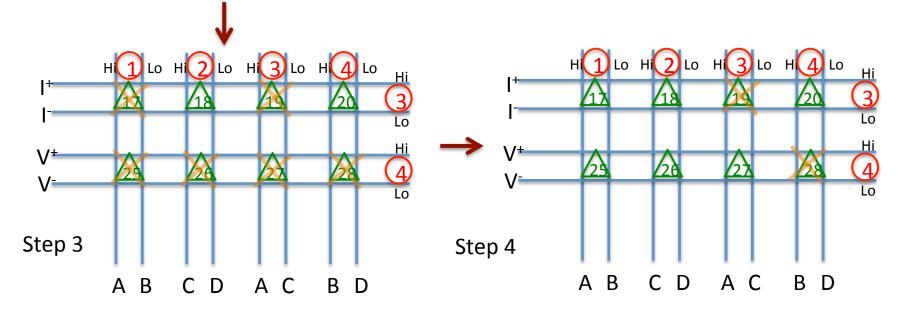


Steps 5, 6, 7 Measurement Sequence – R_{AC, BD}

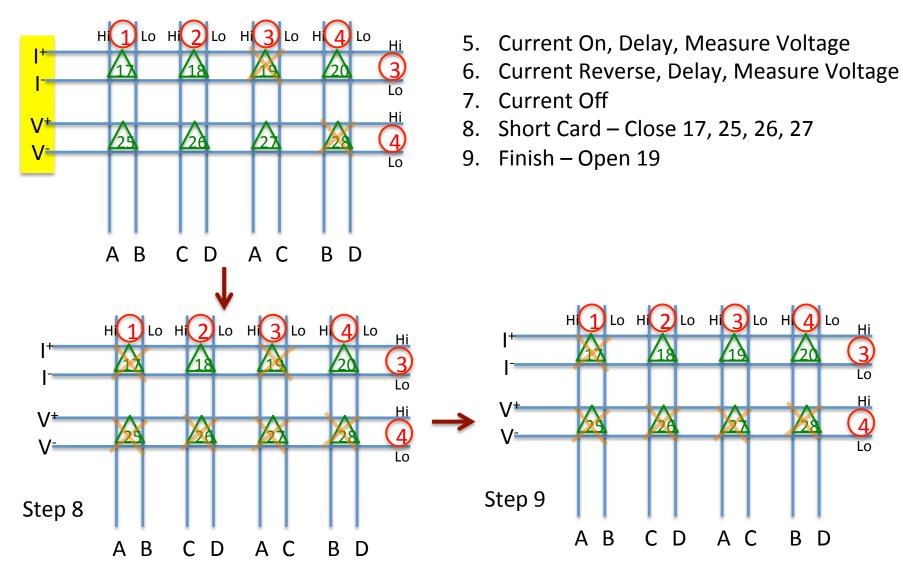

Steps 1, 2 Measurement Sequence – R_{BD, AC}

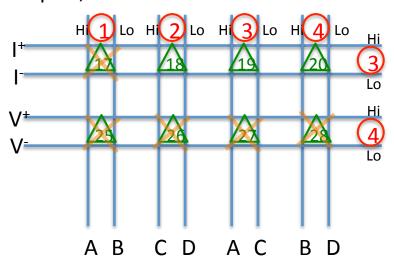

- Current Off
- 2. Short Card Close 17, 25, 26, 27, 28
- 3. Prepare Measurement Close 20
- 4. Prepare Measurement Open 17, 25, 26, 28

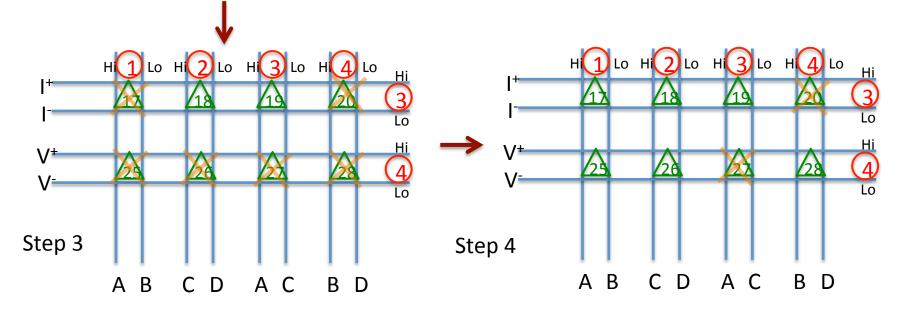

Steps 5, 6, 7 Measurement Sequence – R_{BD, AC}

